创普。 Hailian Packaging Equipment Co.,Ltd mobilemilkingmachine.com

Cow hip lifter, rear hip lift, hip lifting cow bar, cow harness for lifting, cattle hip huggers, lifting device

Basic Information

. Place of Origin: Jiangsu, China(Mainland)

Brand Name: Chuangpu · Certification: ISO HL-Q16 Model Number: Minimum Order Quantity: 1pcs

• Price: USD\$134.00/PCS-USD\$135.00/PCS

CARTON Packaging Details:

• Payment Terms: T/T, L/C, Western Union, MoneyGram

Supply Ability: 1000PCS PER MONTH



Product Specification

Brand Name: Chuangpu Model Number: HL-Q16

Place Of Original: Jiangsu China(Mainland)

Yes Is_customized: · Condition:

Material: Stainless Steel Rubber Galvanized Tube

Tool Net Weight: 8.30kg 600mm Height: • Width: 290mm 700mm · Length:

• Highlight: cattle hoof trimming discs, cattle nose ring



More Images









Product Description

Cow hip lifter, rear hip lift, hip lifting cow bar , cow harness for lifting , cattle hip huggers

Product Model Number: HL-Q16

Product Name: Cow Stand Support, Cow Stand Up Assist Device, Cow Stainless Steel Standing Stand/Holder/Frame


Product Material: Stainless steel, rubber, galvanized tube

Product Technical Data: Net Weight 8.30kg, Length 700mm, Width 290mm, Height 600mm.

Use with a lifting device to lift a cow that is down from milk fever, calving injury or a fall. Rubber-padded stirrups protect hip bones, and movable parts are sealed and lubricated. Adjustable from small heifers to large cows.

1—0···9···	700mm
Width	290mm
Height	600mm
Net Weight	8.30kg



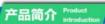

倒普机械 Chuang Pu



Product name: Cow Stand Support

Product number: HL-Q16

Product introduction: The stand support is a kind of machine that cattle can be easily lifted by the use of traction, loading and unloading machines, triangular hoists, etc. when the cattle are fallen, and is suitable for dystocia, severe typhoid fever, fever and other symptoms. The product is made of high-strength iron and steel rods, can carry a weight of 1000KG, and has a rubber tube to protect the cow's buttocks. It is convenient to adjust the spacing by turning the handle.



Installation



First install one end of the column angle



Install the other end of the column



use wire to assist the m to close together and insert bolts to lock



Insert handle




The installation is complete

Product Manual

Place the stand support to the same width as the burdock and place the burdock position. Then shake the stand handle to make it clamp the burdock and lift the cow through the traction equipment. It should be noted that the cow can not be lifted for more than 3 hours to avoid the congestion of the burd bone.



创普机械 Chuang Pu

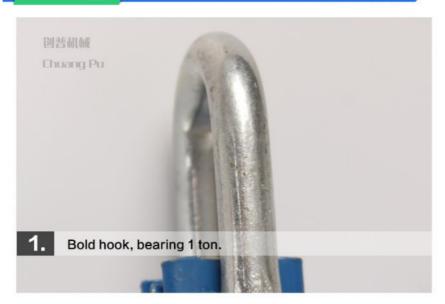












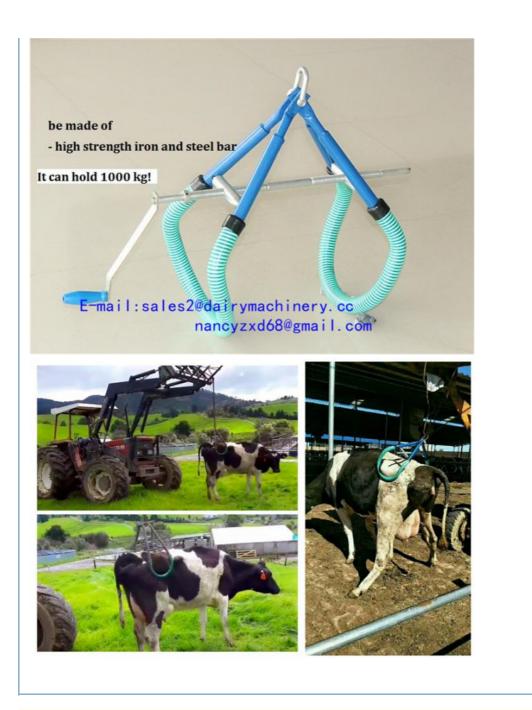



Hailian Packaging Equipment Co.,Ltd



0086 13301569240



sales2@dairymachinery.cc



@ mobilemilkingmachine.com

Nanzhuang Village, Economic Development Zone, Yangshe Twon, Zhangjiagang City, China