

## Steel Cow Horn Cutter Machine Dehorner For Dairy Farm With 1.7KW LEIYA Motor

Our Product Introduction

for more products please visit us on [mobilemilkingmachine.com](http://mobilemilkingmachine.com)


### Basic Information

- Place of Origin: jiangsu
- Brand Name: chuangu
- Certification: ISO CE FDA SGS
- Model Number: HL-Q15A
- Minimum Order Quantity: 5 sets
- Price: USD\$680-690/set
- Packaging Details: carton
- Delivery Time: 2-4days
- Payment Terms: T/T, Western Union, Paypal
- Supply Ability: 500sets/month



### Product Specification

- Voltage: 220V/50hz
- Power: 1.7KW
- Package: Carton
- Usage: Horn Cutting
- Rotating Speed: 8500r/Min
- Weight: 8.6KG
- Highlight: cattle hoof trimming discs, cattle nose ring

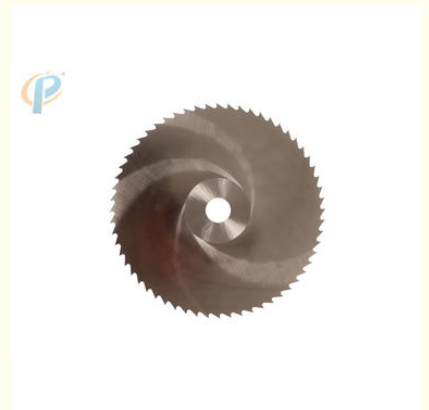

### More Images



## Product Description

### Product introduction:

Cow Dehorner, cutting the horn annular portion from the its bottom with an electric cutter after holding cows, fast dehorning and easy operation. Removing horn prevents cow from getting hurt in fights, especially the breast part is not gored by the following cow horn. The cows after dehorning are quieter and easier to manage.

### Product information

Model: HL-Q15A

Item: Horn Cutter

Usage: Horn cutting for cow

Voltage	220v/50hz
Power	1.7KW
Package	Carton
Rotating Speed	8500r/min
Power Cord Length	2.25m
MOQ	5 sets

## Product parameters

创普机械  
Chuang Pu



型 号 Model	电压 Voltage	功率 Power	材质 Material
HL-Q15A	220V/50HZ	1700W	钢 ( steel )
转速 Rotating speed	刀盘直径 Plate diameter	电源线长 Power cord length	重量 Weight
8500RPM/MIN	225MM	2.25M	8.6KG



Weight: 11KG

# How to use the Cow Dehorner



1. Fix cow head on the cow headlock other devices with fixed functions with a tool such as a rope.
2. Turn on the power, and or place the annular cutting hole of the dehorner at the horn on the bull's head.
3. By pushing the auxiliary handle, the horn can be cutted instantly by saw blade and cow hardly feels any pain.

## Attention!

As to young cow (not pregnant or in second trimester) and dairy cows in the late milking period, to avoid causing miscarriage and a sharp drop in milk production, the dehorning time should be selected in the seasons with less mosquitoes and cool weather to prevent infection.

During the process of cutting horn, the saw blade rotates at high speed and produces extremely high temperature, and then the capillaries in the horn have been burnt, so this method has less bleeding and faster speed. If there are bleeding capillaries, it can use an electric iron to stop bleeding, and then anti-inflammatory drugs are used to disinfect the wound.

The horn cavity is filled with sterilized cotton after removing horn remember not to plug it too tightly. Generally, the wound can be naturally flat for about two months. Please observe the condition of cow after dehorning, and deal with the infection in time.



**Hailian Packaging Equipment Co., Ltd**



0086 13301569240



sales2@dairymachinery.cc



mobilemilkingmachine.com

Nanzhuang Village, Economic Development Zone, Yangshe Twon, Zhangjiagang City, China